

Matrix® Hardware

Matrix Easy Set hardware is easily installed in minutes and all adjustments can be made with the user seated in the wheelchair.



Features and Functional Benefits

- Quick release for easy portability
- Infinite adjustment of depth - (4") angle and 40° rotation
- Fits 3/4", 7/8" and 1" diameter canes

Hardware Adjustments



DEPTH:
2" behind back canes



DEPTH:
2" in front of back canes



ANGLE:
20° posterior range
of back angle adjustment



ANGLE:
20° anterior range
of back angle adjustment



ROTATION:
Back may be oriented to
accommodate pelvic/trunk rotation

HEAVY DUTY "EASY SET" MOUNTING HARDWARE



Standard with Heavy Duty Posture Backs.
Recommended for users over 130kg, high tone or heavy use.

FIXED SPACE SAVER MOUNTING HARDWARE




Requires tools for removal.
Saves space for tighter width applications
Fully adjustable.

SPLIT MOUNTING HARDWARE



Provides clearance of mid back cane rigidiser bars.
Not recommended for Heavy Duty backs.

EXTENDED MOUNTING PINS

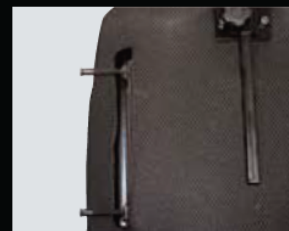

Allows for back to fit chairs up To 2" wider.
Not available for Heavy Duty Backs.



0800 212 333 email: info@c1south.co.nz
www.c1south.co.nz